

LINSEIS


T H E R M A L A N A L Y S I S

DILATOMETRY

DIL L76

DIL L75 Horizontal

DIL L75 Vertical

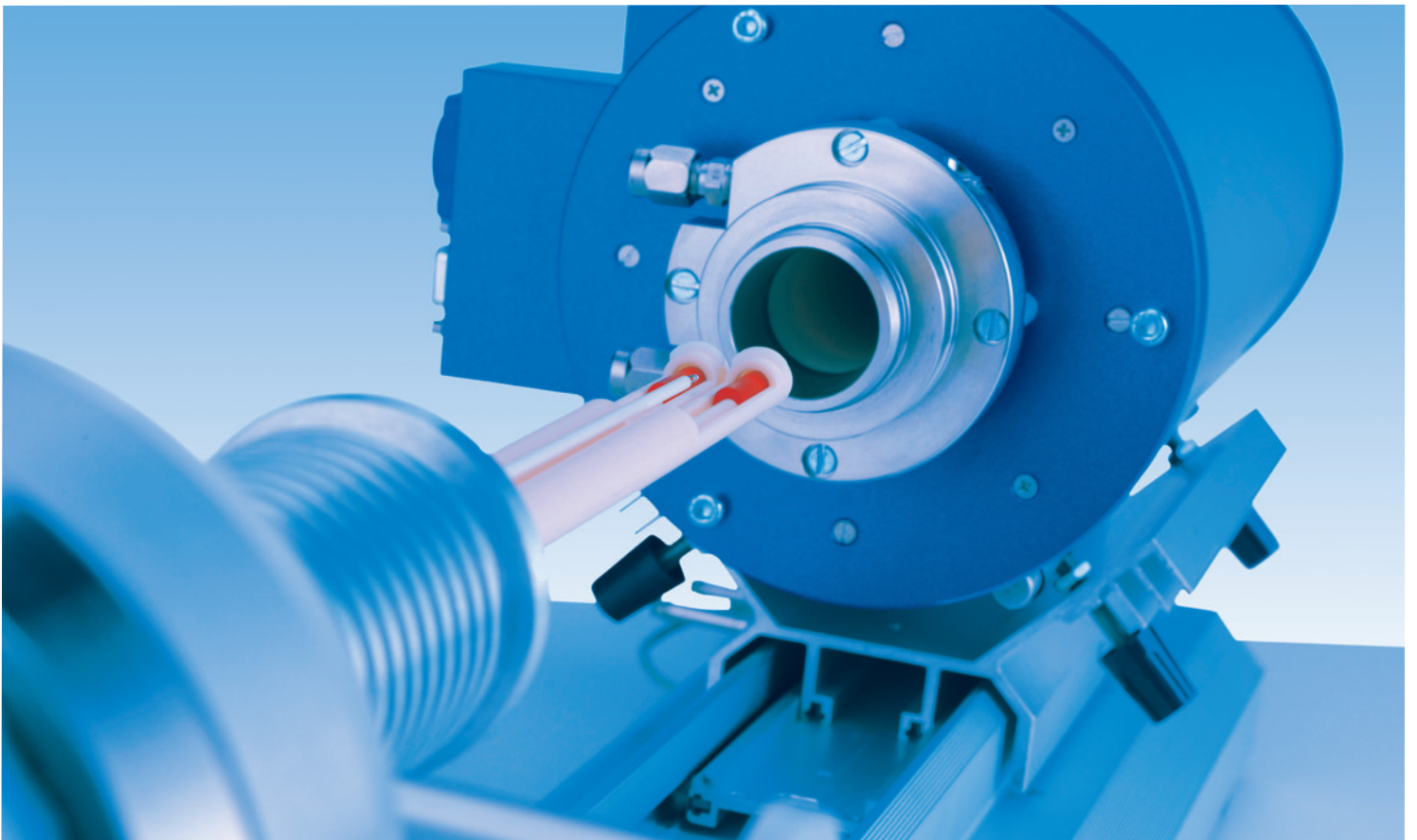

Since 1957 LINSEIS company has been delivering outstanding service, know how and leading innovative products in the field of thermal analysis and thermo physical properties.

Customer satisfaction, innovation, flexibility and high quality are what LINSEIS represents. Thanks to these fundamentals our company enjoys an exceptional reputation among the leading scientific and industrial organizations. LINSEIS has been offering highly innovative benchmark products for many years.


The LINSEIS business unit of thermal analysis is involved in the complete range of thermo analytical equipment for R&D as well as quality control. We support applications in sectors such as polymers, chemical industry, inorganic building materials and environmental analytics. In addition, thermo physical properties of solids, liquids and melts can be analyzed.

LINSEIS provides technological leadership. We develop and manufacture thermo analytic and thermo physical testing equipment to the highest standards and precision. Due to our innovative drive and precision, we are a leading manufacturer of thermal Analysis equipment.

The development of thermo analytical testing machines requires significant research and a high degree of precision. LINSEIS Corp. invests in this research to the benefit of our customers.




Claus Linseis
Managing Director

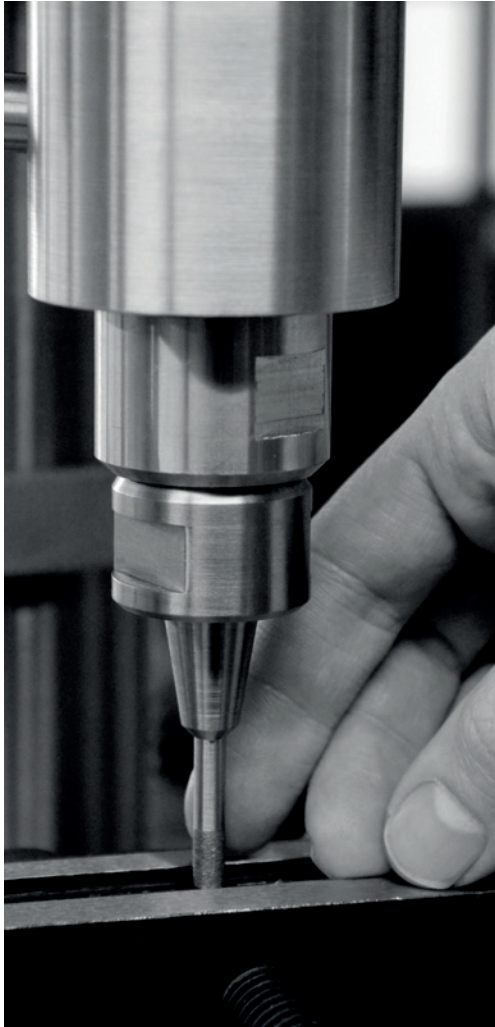

German engineering

The strive for the best due diligence and accountability is part of our DNA. Our history is affected by German engineering and strict quality control.

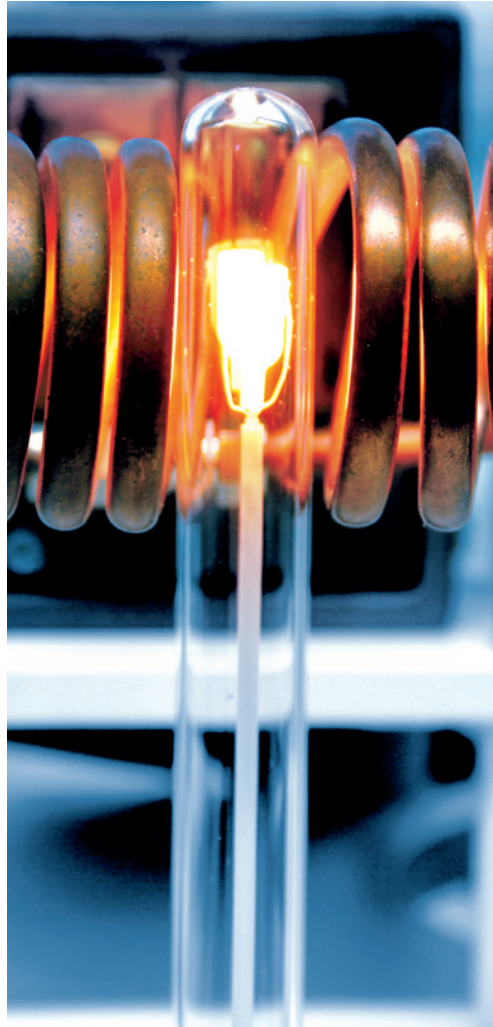

Innovation

We want to deliver the latest and best technology for our customers. LINSEIS continues to innovate and enhance our existing thermal analyzers. Our goal is constantly develop new technologies to enable continued discovery in science.

GENERAL


Dilatometry (DIL) is a technique in which a dimension of a substance (such as: ceramics, glasses, metals, composites, polymers and others) under no or negligible load is measured (e.g. expansion measurement or shrinkage measurement) as a function of temperature while the substance is subjected to a controlled temperature program in a specified atmosphere.

LINSEIS produces a wide range of Single-, Dual-, Differential-, Quattro-, Quenching-, Laser- and Optical Dilatometers in a temperature range from -263°C up to 2800°C and a resolution up to 0.03 nm.

Thanks to our years of experience (est. 1957), Linseis offers Dilatometers with an unbeaten precision variability and performance.

Dilatometers are frequently used for R&D and quality control of solids, liquids, powders and pastes to determine their:

- determination of thermal expansion coefficient (CTE)
- linear thermal expansion (ΔL)
- sinter-temperatures and sinter steps
- determination of glass transition (T_g)
- phase changes
- optimization of burning processes
- volume changes
- Rate Controlled Sintering (RCS)
- decomposition
- density change

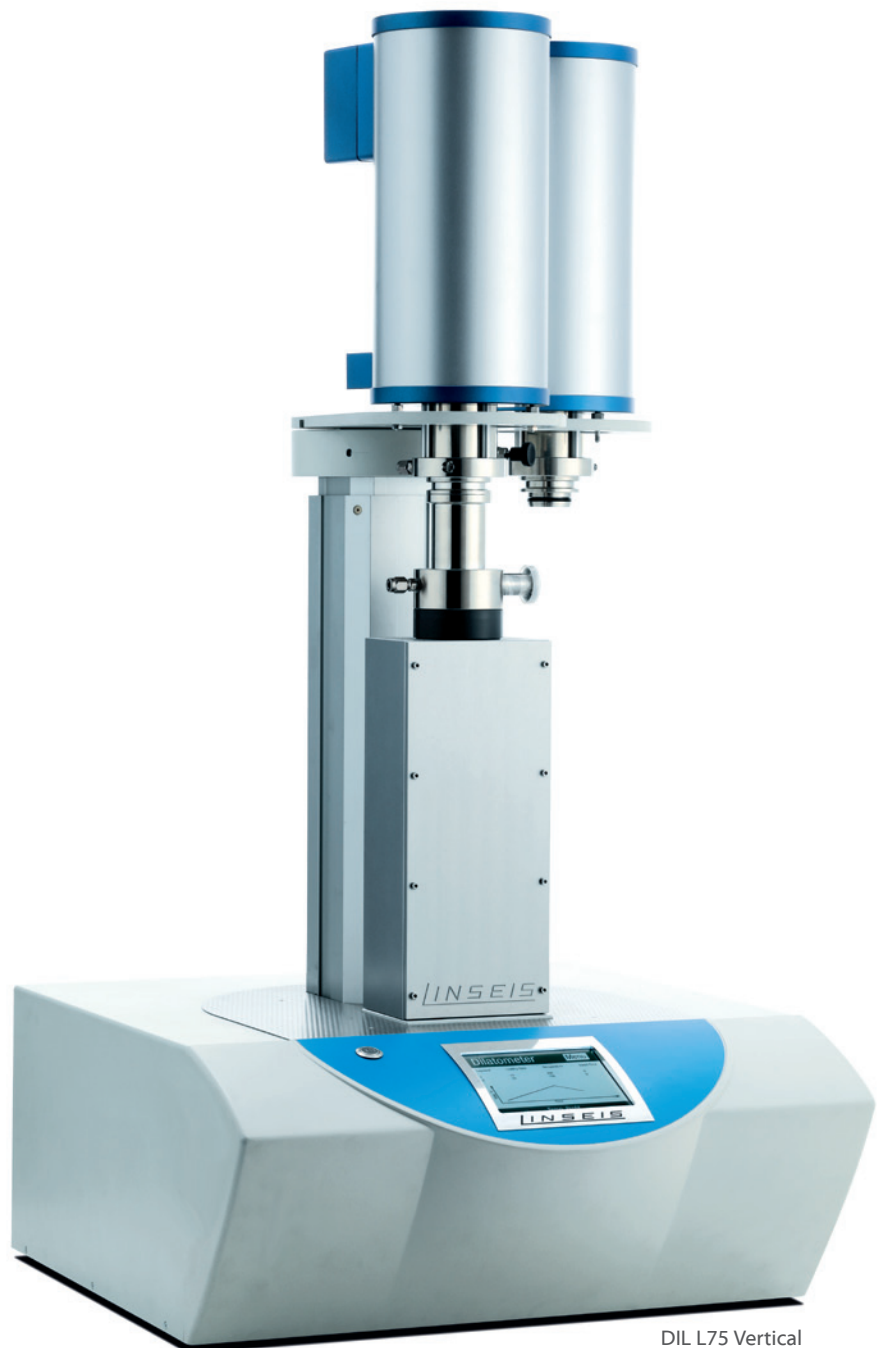

DIL L75 Vertical

Depending on your application, we offer the perfect solution:

Horizontal Dilatometer L76/L75Horizontal:

- multipurpose system
- highest temperature uniformity
- L75H is perfect for research & development

Vertical Dilatometer L75Vertical:

- friction free sample holder
- push-rod contact is always guaranteed
- possible field of application Rate Controlled Sintering (RCS)
- best arrangement for low -263°C and high temperature $>2800^{\circ}\text{C}$ applications


Optical Dilatometer L74 Optical*:

- non-contact design
- only choice for melting point determination
- best Choice for Rate Controlled Sintering
- can be used for additional applications:
contact angle determination and heating microscopy
- best choice for irregularly shaped and soft samples

Quenching and Deformation Dilatometer

L78 Q/D/T*:

- unreachd heating and cooling rates up to 4000K/s
- creation of CHT/CCT/TTT diagrams
- deformation up to 25 kN and tension forces upto 23 kN

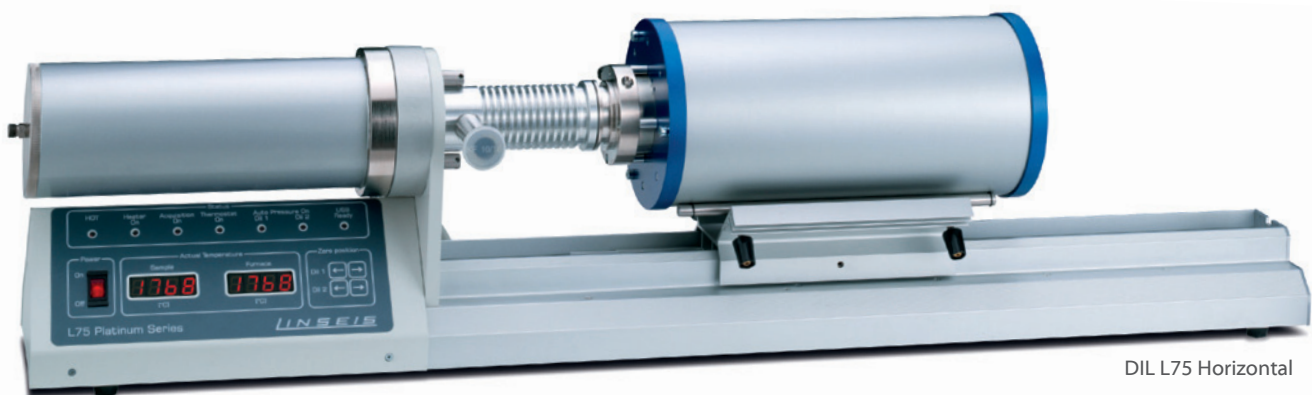

DIL L75 Horizontal


*see separate product brochure



HORIZONTAL VERSUS VERTICAL DILATOMETER

Horizontal arrangement

- easy and robust design at a modest price
- from -263 to 2800°C. For applications requiring lower or higher temperatures vertical arrangement is better suited

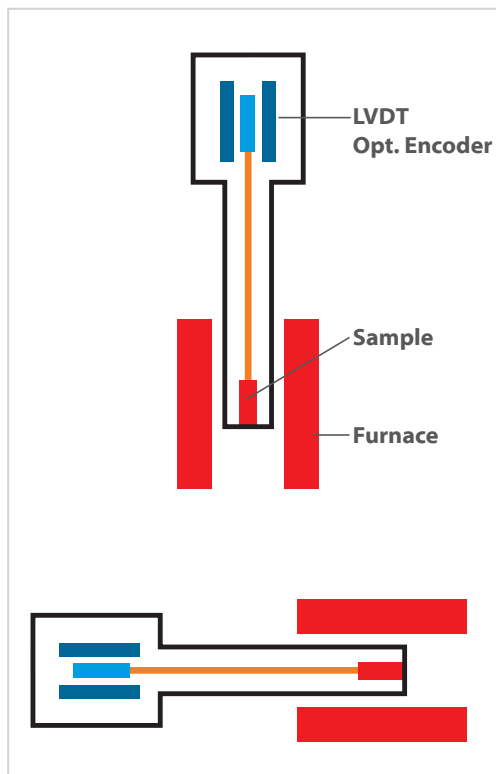

Vertical arrangement

- zero friction design (sample is only in contact with the end knob and the pushrod)
- multi-furnace arrangement (up to three furnaces)...
- to cover entire temperature range from -180 to 2800°C
- to increase throughput (hot furnace can be lifted automatically to switch to new cold furnace and start new measurement)
- small footprint
- available as single-, differential/double- and quattro arrangement (1, 2 or 4 samples at the same time)
- perfect for low temperature measurements (furnace on the bottom - measurement compartment on top) to ensure natural gas path, cold air (falls down) below sensor compartment
- perfect for high temperature measurements (furnace on top - measurement compartment below) to ensure natural gas path, hot air (streams up) above sensor compartment

COMPARING LVDT AND LINEAR OPTICAL ENCODER

LVDT

The LVDT (Linear Variable Differential Transformer) consists out of 3 coils: the LVDT body and the movable core. The primary coil is excited by a low frequency AC voltage. The two secondary coils are put in series with reversed polarity. The horizontal position of the core defines the amount of coupling between the primary and the secondary coils. When the core is in the center position, the induced voltages in the secondary coils have the same amplitude. Due


to reversed polarity of the two coils the sum (output voltage) is zero. When the core moves, the coupling between primary and secondary coils changes. So, in one secondary coil the induced voltages increases, while in the other the voltage decreases. So, the sum of both is no longer zero. The sum of amplitude depends on the amount of movement of the core, while the phase (polarity) depends on the direction of the movement.

Advantage:

- the output signal is absolute and unique for each position, no reference movement necessary
- the core can be moved without any friction
- the resolution is infinitive, limited by the noise of the electronics used for signal conditioning
- not sensitive better suited for applications in dirty environments (gas, vacuum, dust)

Disadvantage:

- limited measuring range, e.g. +/- 2,5mm
- needs calibration

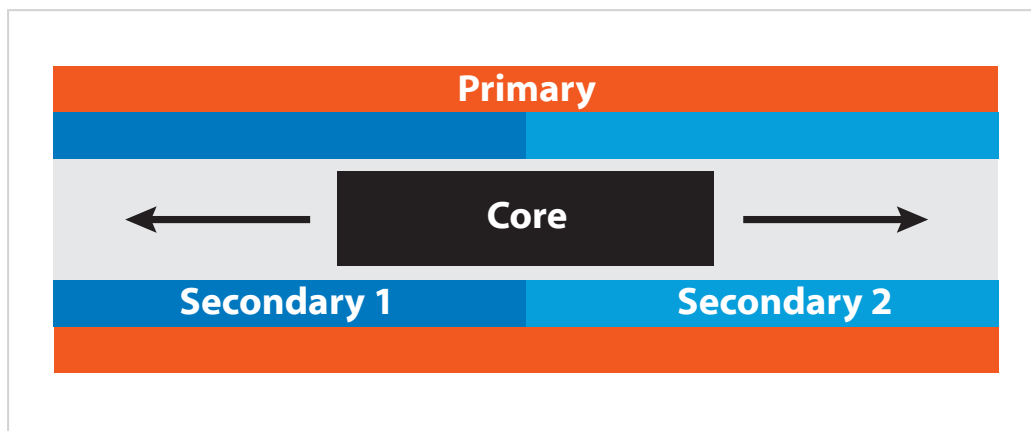

Linear Optical Encoder

A linear optical encoder uses a ruler, made of glass or metal, with a special optical pattern on it. Usually transparent and non-transparent, or reflecting and non-reflecting lines are used. A light source shines to the ruler, and the bright/


dark transitions are measured. The count of measured transitions corresponds to the displacement, while the phase contribution of the two detectors A and B depends on the direction of the movement.

Advantage:

- no calibration necessary
- no limit of measuring range
- force modulation – offering TMA capabilities
- superior reproducibility

Disadvantage:

- the position change is measured relative, a zero position measurement for absolute position reading is necessary
- limited resolution: the lower limit for the distance of the patterns on the ruler is app. 20nm. For higher resolutions the reading must be interpolated
- the optical detector system is sensitive to dust (production environment)



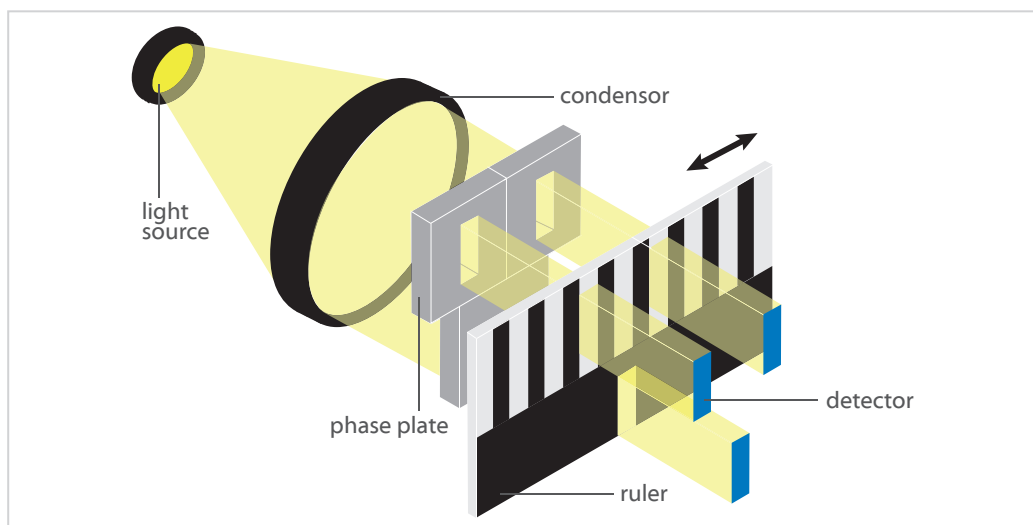


FEATURES

Exciting new features to simplify work.

Measuring systems

Numerous different single, double or quattro measuring systems made of fused silica, Al_2O_3 or graphite are available.

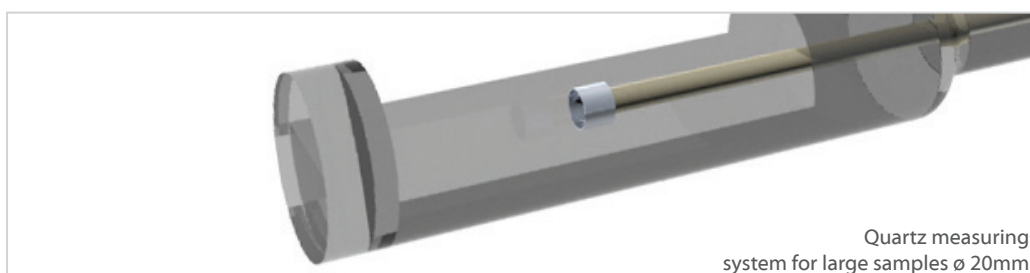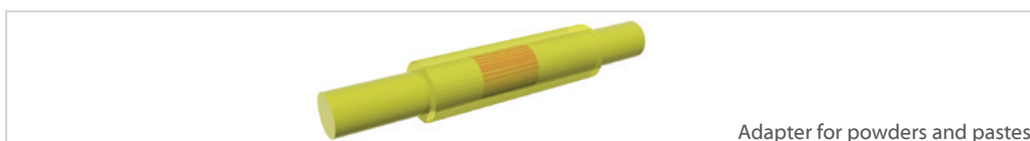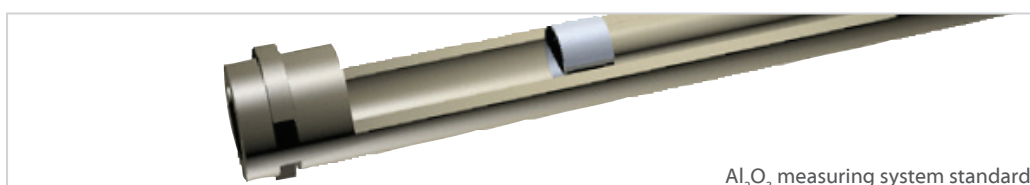

Instrument design

Optimum sensor selection to ensure best measurement result for any given requirement. Integrates Mass Flow Controllers or gas dosing system. Thermostatically controlled measuring compartment disables environmental measurement influences.


MultiNudge

Our sample positioning feature ensures an optimum sample alignment. As the sample piston touches the sample several times it improves

the orientation of the sample which is utmost important for a successful dilatometer measurement.

Automatic Contact Force Adjustment

When it comes to samples with large expansion a contact force adjustment during the measurement is critical to ensure excellent measuring results. This feature eliminates the influence of varying force during the measurement.

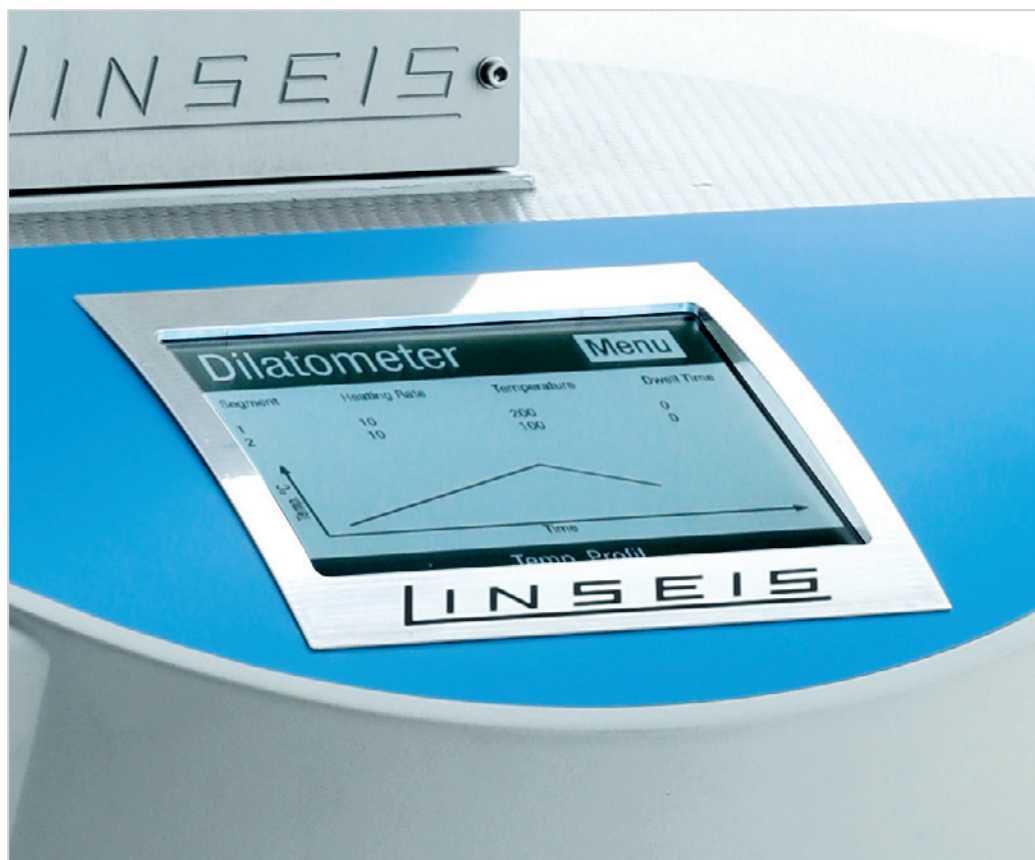

Vacuum and controlled atmosphere

The instrument design provide for high vacuum, inert, reducing, oxidizing or humidified atmosphere. Furthermore, the instrument can be pressurized up to 5 bar overpressure (option). Measurements can be done under corrosive conditions if proper conditions are respected. Residual gas analysis systems can be coupled via an optional heated capillary

Wide temperature range -263 to 2800°C

The LINSEIS Dilatometer instruments can be equipped with up to three furnaces at the same time. A broad variety of different furnaces is available to enable measurements in the widest temperature range on the market. Unmatched selection of furnaces for widest possible temperature range.

Evolved gas analysis

Optional gas analysis with MS, FTIR or GCMS is possible. This provides valuable additional information of the sample composition.



SOFTWARE

All LINSEIS thermo analytical instruments are PC controlled. The individual software modules exclusively run under Microsoft® Windows® operating systems.

The complete software consists of 3 modules: temperature control, data acquisition and data evaluation. The software incorporates all essential features for measurement preparation, execution, and evaluation of a dilatometer run.

Thanks to our specialists and application experts, LINSEIS offers comprehensive easy to understand user friendly application software.

Software-Features:

- program capable of text editing
- data security in case of power failure
- thermocouple break detection
- repetition measurements with minimum parameter input
- evaluation of current measurement
- curve comparison


Software

- storage and export of evaluations
- export and import of data ASCII
- data export to MS Excel
- multi-methods analysis (DSC TG, TMA, DIL, etc.)
- zoom function
- 1st and 2nd derivative

- programmable gas control
- statistical evaluation package
- automatic axis re-scaling
- softening point detection

DIL Features:


- Rate Controlled Sintering (RCS) software
- interchangeable thermocouples for various atmospheres
- sinter process evaluation
- glass transition and softening point evaluation
- softening point determination and system shut down
- linear thermal expansion evaluation
- several system correction features
- automatic zero point adjustment
- auto-scheduler for up to 16 sub-sequent runs

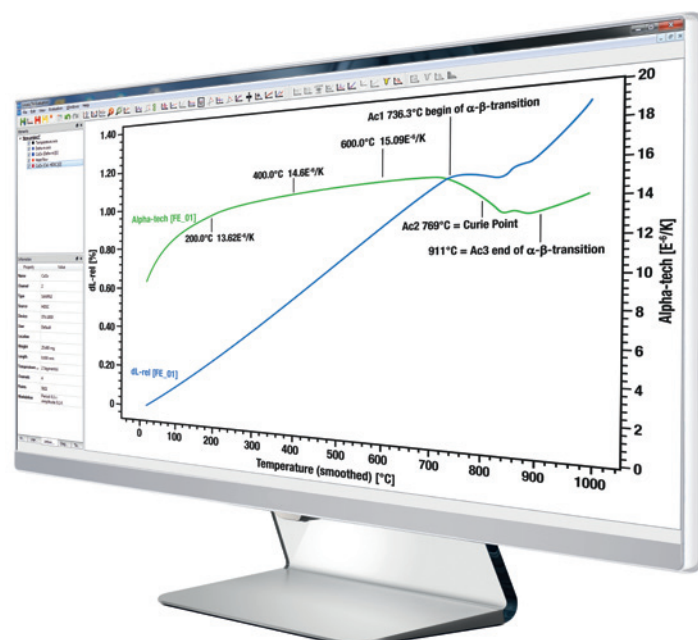

Thermal Library

The LINSEIS thermal library software package comes as an option for the well-known, user friendly LINSEIS Platinum evaluation software, that is integrated in almost all our instruments.

The thermal library allows for the comparison of the complete curves with a data base providing thousands of references and standard materials within only 1-2 seconds.

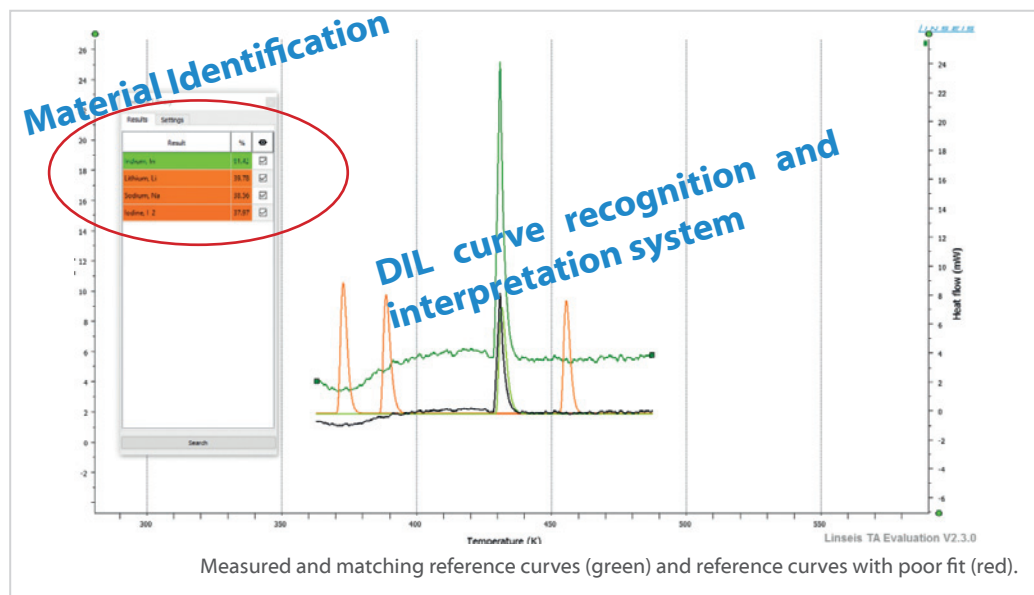

Multi-Instrument

All LINSEIS instruments DSC, DIL, STA, HFM, LFA, etc. can be controlled from one software template.

Multi-User

The administrator can generate different user levels providing different rights to operate the instrument. An optional log file is available, too.

Report Generator

Convenient template selection to generate customized measurement reports.

Kinetic software

Kinetic analysis of DSC, DTA, TGA data for the study of the thermal behavior of raw materials and products.

Multi-Lingual

Our software is available in many different user exchangeable languages, such as: English, Spanish, French, German, Chinese, Korean, Japanese, etc.

Data Base

Physical constants of new materials can be easily integrated by the user.


SPECIFICATIONS

	L76	L75 Horizontal	L75 Vertical	
Temperature range	RT up to 1600°C	-180 up to 2800°C	-263 up to 2800°C	
Delta L resolution	0.05 nm	0.03 nm	0.03 nm	LVDT
Measuring range	±2500 µm	±2500 µm	±2500 µm	
Contact force	—	10 mN up to 1N	10 mN up to 1N	
Delta L resolution	1 nm	0.1 nm	0.1 nm	Optical Encoder
Measuring range	±25000 µm	±25000 µm	±25000 µm	
Automatic sample length detection	yes	yes	yes	
Force modulation	no	yes	yes	
Contact force	50mN up to 3N	10 mN up to 5N	10 mN up to 5N	
Heating rates	Based on furnace: <ul style="list-style-type: none"> • steel, copper, fused silica, silicon carbide: 0.001 up to 50K/min • graphite 0.001 up to 100 K/min 			
Sample holder	user interchangeable, SiO ₂ , Al ₂ O ₃ , graphite All sample holders available as: <ul style="list-style-type: none"> • single sytem (one pushrod) • system with two pushrods usable in dual or differential mode • Al₂O₃ sample holder SiO ₂ and Al ₂ O ₃ sample holders available as hollow tubes or massive rods			
Temperature /Accuracy / Precision / Resolution	1 K / 0.1 K / 0.001 K	1 K / 0.1 K / 0.001 K	1 K / 0.1 K / 0.001 K	
Thermal stability (isotherm)	± 0.02 K	± 0.02 K	± 0.02 K	
ΔL / L₀ Repeatability	0.001%	0.001%	0.001%	
ΔL / L₀ Accuracy	0.002%	0.002%	0.002%	
Force resolution	—	0.001 mN	0.001 mN	
Gas atmosphere	Inert, oxid., red.	Inert, oxid., red., vac.	Inert, oxid., red., vac.	
Software	latest Windows operating system older versions possible on request			

	L76	L75 Horizontal	L75 Vertical
Multiple furnace configuration	optional	up to 2 furnaces	up to 3 furnaces
Motorized furnace operation	optional	optional	included
Gas dosing	manual gas dosing or mass flow controller 1/3 or more gases		
Contact force adjustment	—	included	included
Single / double dilatometer	—	optional	optional
Softening point detection	included	included	included
Density determination	optional	included	included
L-DTA	optional	included (up to 2000°C)	included (up to 2000°C)
Rate Controlled Sintering (RCS)	optional	included	included
Thermal library	optional	included	included
Electric thermostization of measuring head	included	included	included
Low temperature options	LN ₂	LN ₂ , Intra	LN ₂ , Vortex, Helium cryostat
Vacuum tight design	—	yes	yes
Automatic evacuation system	—	optional	optional
OGS Oxygen Getter System	—	optional	optional


FURNACES

Temperature	Type	Heating element	Atmosphere	Temperature sensor
-263 – 300°C	L75/264 He	Thermo coax	inert, oxid., red., vac.	Semiconductor / Pt 100
-180 – 500°C	L75/264	Thermo coax	inert, oxid., red., vac.	Type K
-180 – 700°C	L75/264/700	Thermo coax	inert, oxid., red., vac.	Type K
-180 up to 1000	L75/264/1000	Thermo coax	inert, oxid., red., vac.	Type K
RT – 1000°C	L75/220	Kanthal	inert, oxid., red., vac.	Type K
RT – 1400°C	L75/230	Kanthal	inert, oxid., red., vac.	Type S
RT – 1600°C	L75/240	SiC	inert, oxid., red., vac.	Type S
RT – 1650°C	L75/240 PT	Platinum	inert, oxid., red., vac.	Type S
RT – 1750°C	L75/240 M	MoSi ₂	inert, oxid., red., vac.	Type B
RT – 2000°C	L75/260	Graphite	N ₂ /Vac.	Type C and/or pyrometer
RT – 2400°C	L75/270	Graphite	N ₂ /Vac.	Pyrometer
RT – 2800°C	L75/280	Graphite	N ₂ /Vac.	Pyrometer

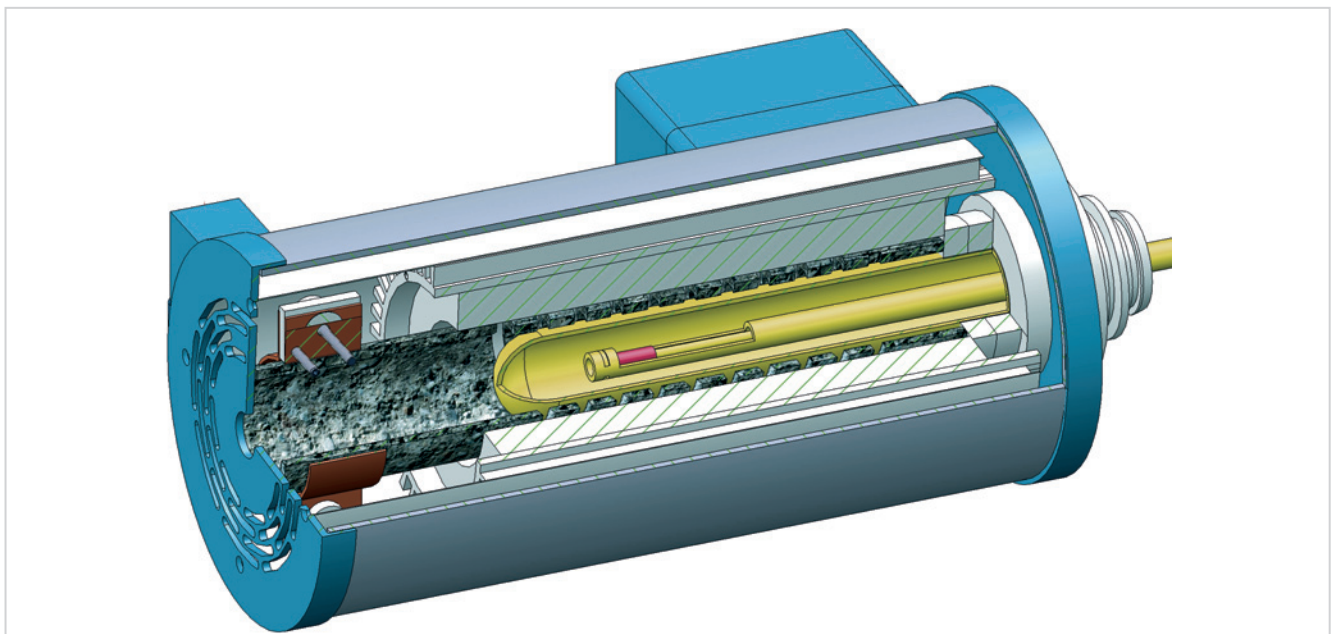

APPLICATIONS


Glass Ceramic

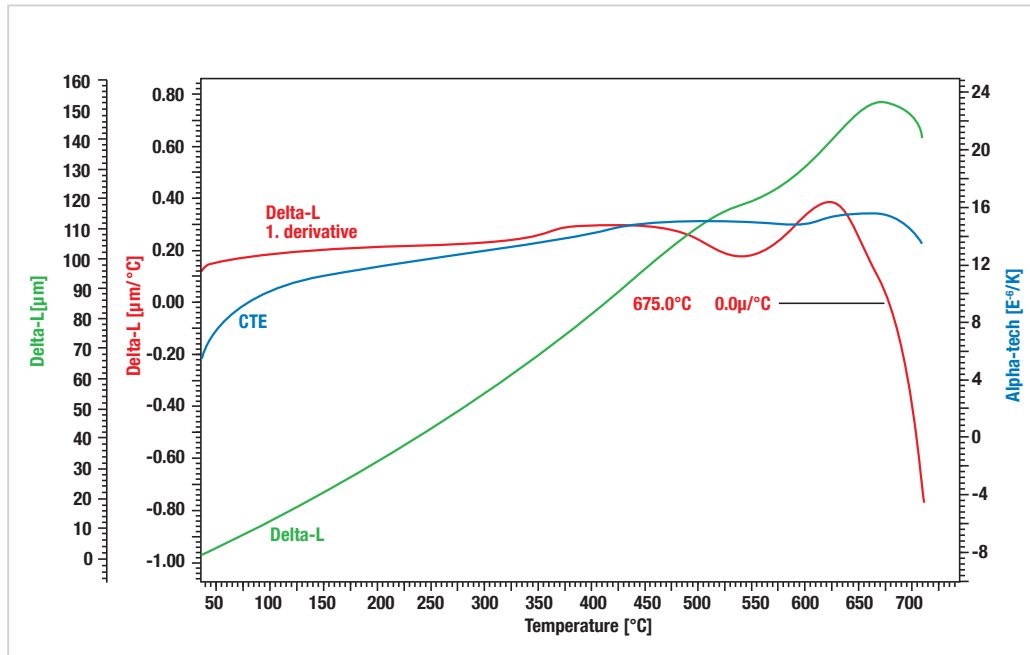

The dilatometric method is an excellent method to determine the thermal expansion (CTE) and the softening point of glass ceramic materials. Beside the absolute expansion and the expansion coefficient (CTE), you can find the first derivative of the absolute expansion. The temperature at which the derivative becomes zero corresponds to the maximum length of the sample and to the softening point of the material.




DTA - Feature

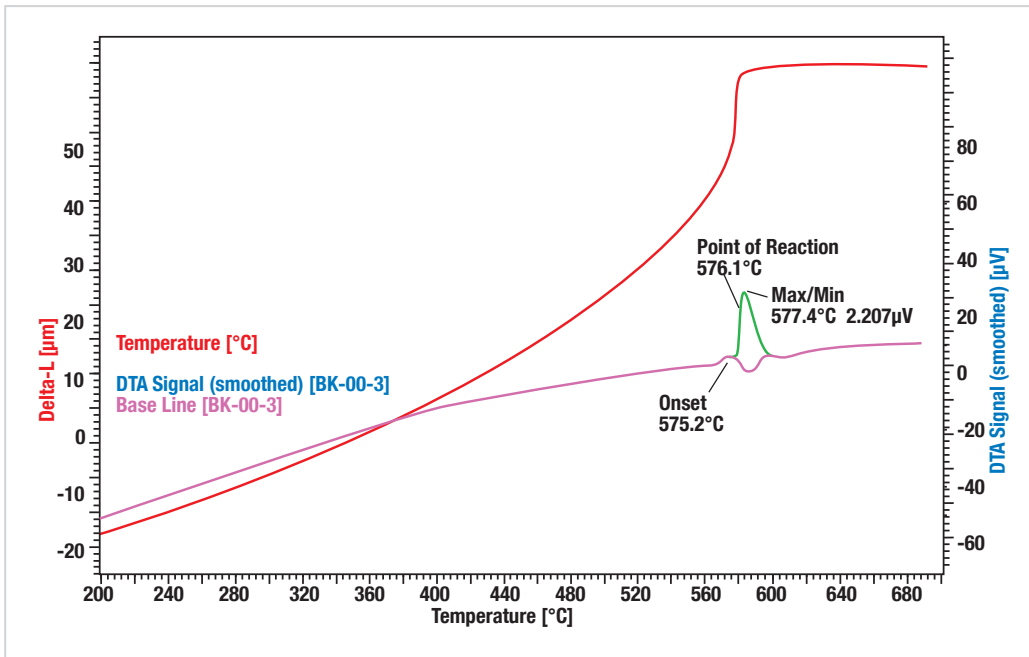


The thermal expansion of α -SiO₂ can be easily evaluated with the L75 Dilatometer. The additional DTA feature enables an in depth view of the thermal behavior of the material. The DTA measurement is a mathematical routine based on the sample temperature. Exo- and endothermic effects influence the change of the sample temperature during the dynamic heating or cooling cycle. The phase transition from α - to β -SiO₂ takes place at app. 575°C. The deviation of the measured temperature from the literature value (574°C) can be used for a temperature calibration.

Glass Sample, softening point

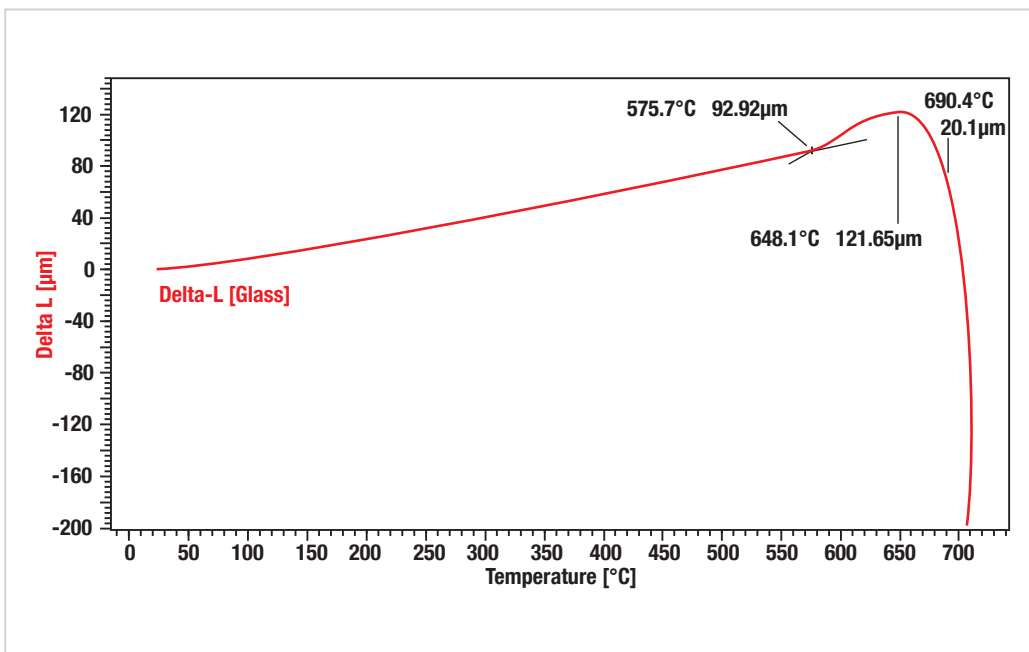

The glass sample was measured using the softening-detection feature. The heating step will be stopped by the software either if the target temperature is reached or if a shrinkage above a pre-defined threshold is detected. In this example, the threshold was set to -100µm meaning that in case that the sample length falls by more than 100µm below the maximum length of the sample, the heating was interrupted and a cooling segment was started automatically. The softening point can be detected in an easy and safe procedure. The automatic stop of the heating avoids damages on the instruments due to sticking of the sample on the measuring system.


Iron

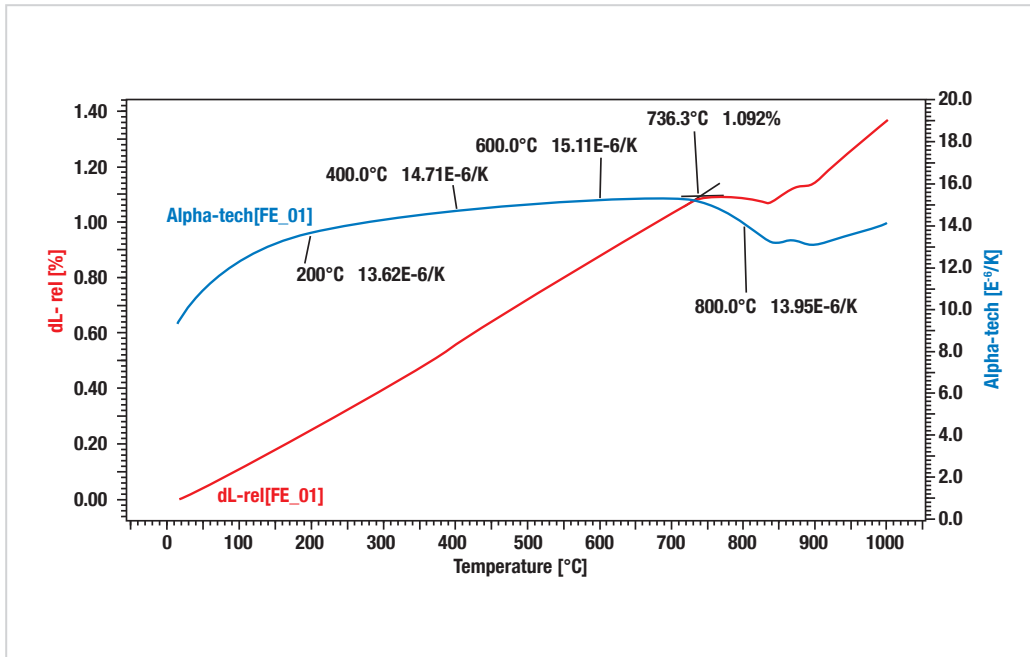

The linear thermal expansion (ΔL) and the CTE of the iron sample under argon atmosphere are evaluated. The heating rate was 5K/min. Above 736.3°C (maximum temperature of CTE) shrinkage was detected, which is due to a change in the atomic structure, known as the Curie-point. The difference of measured and literature value can be attributed to contamination of the sample.




Ceramics / Powder metallurgy

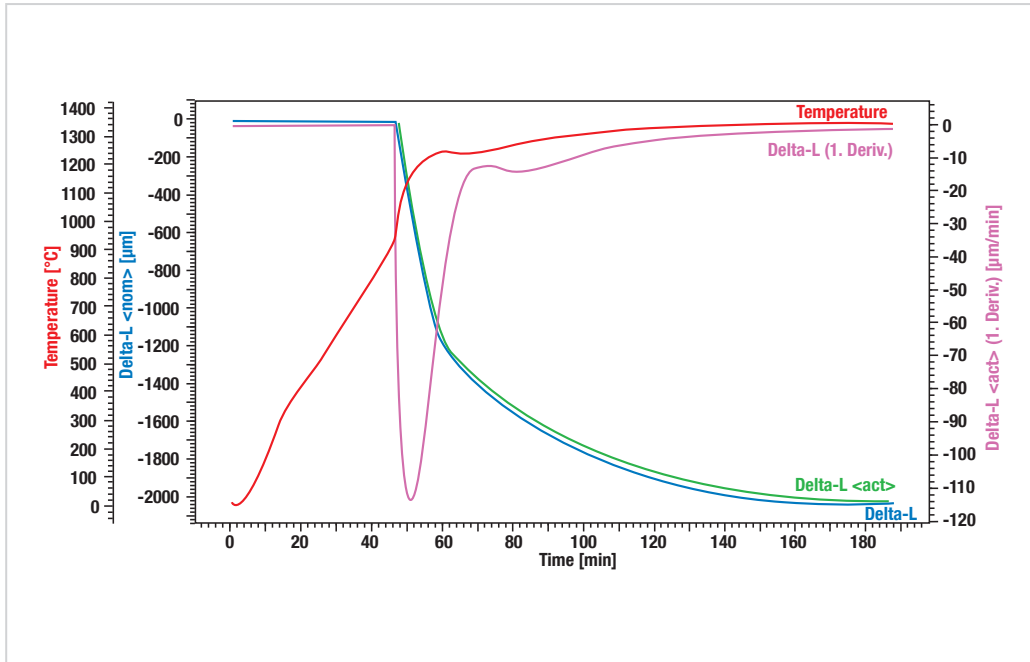

In production processes of high-tech ceramics a simulation of sinter processes is of high interest. When using the optional software package RCS (Rate Controlled Sintering) it is possible to program controlled sintering with a dilatometer as per PALMOUR III theory. The following application is a sintering process of ZrO_2 . After approximately 60 minutes the main sintering step is finished and the heating rate can be slowed down. After approximately 160 minutes the final density of the material is reached and sintering finished.


Copper

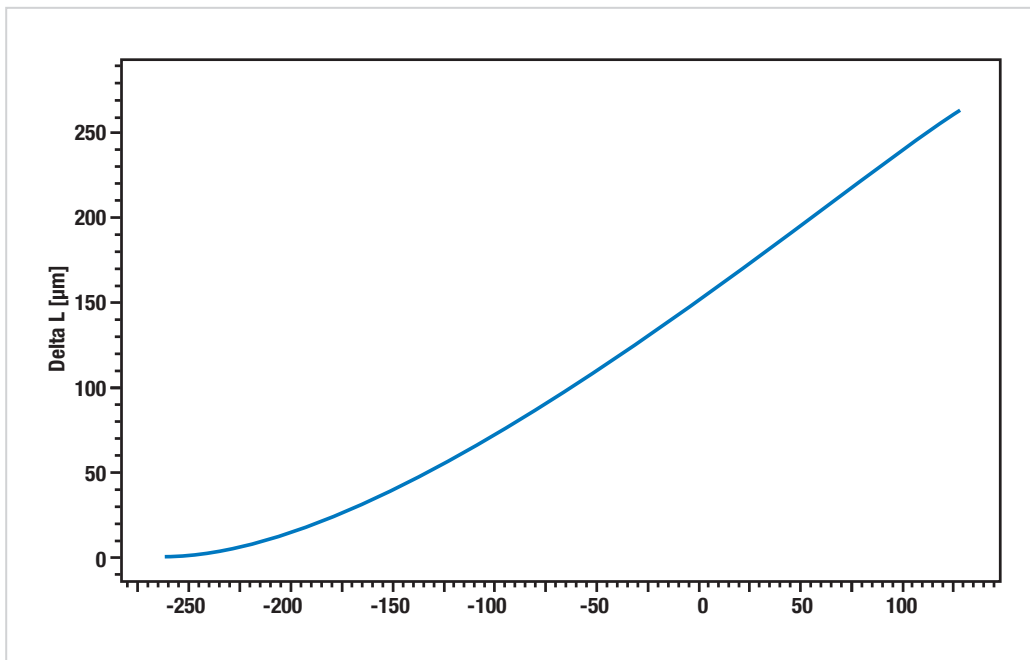

Expansion measurement of copper down to a temperature of 10K using a Helium cryostat.


ALTERNATIVE DILATOMETERS

Optical versus Pushrod-Dilatometer

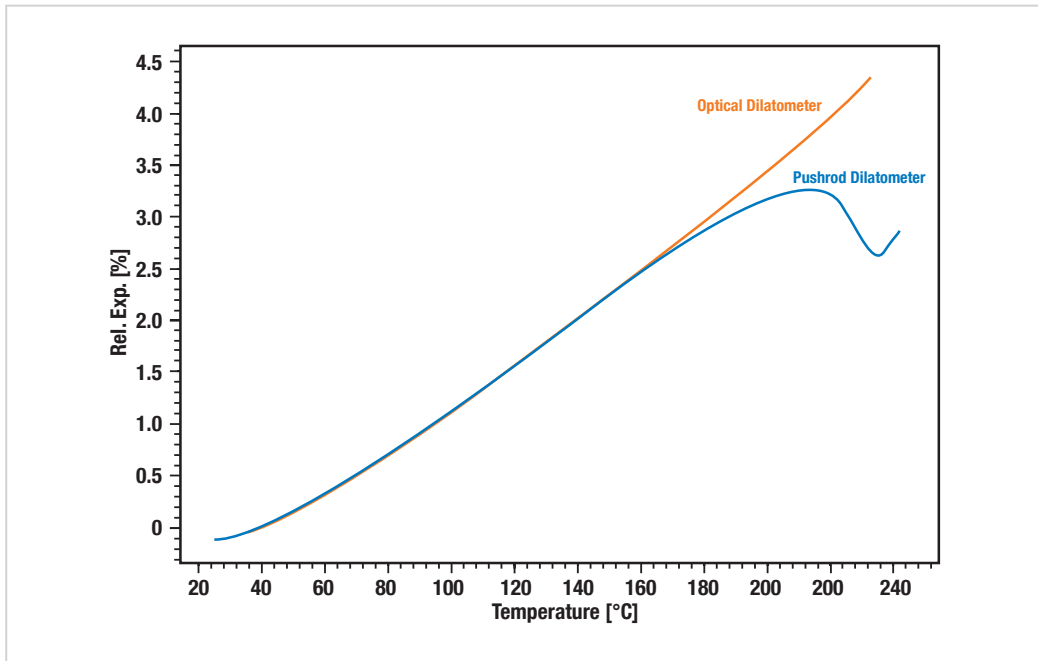

Comparison between Conventional and Optical Dilatometer when evaluating the expansion of an epoxy resin into the melting stage.

(please look at our separate L74 Optical DIL brochure)

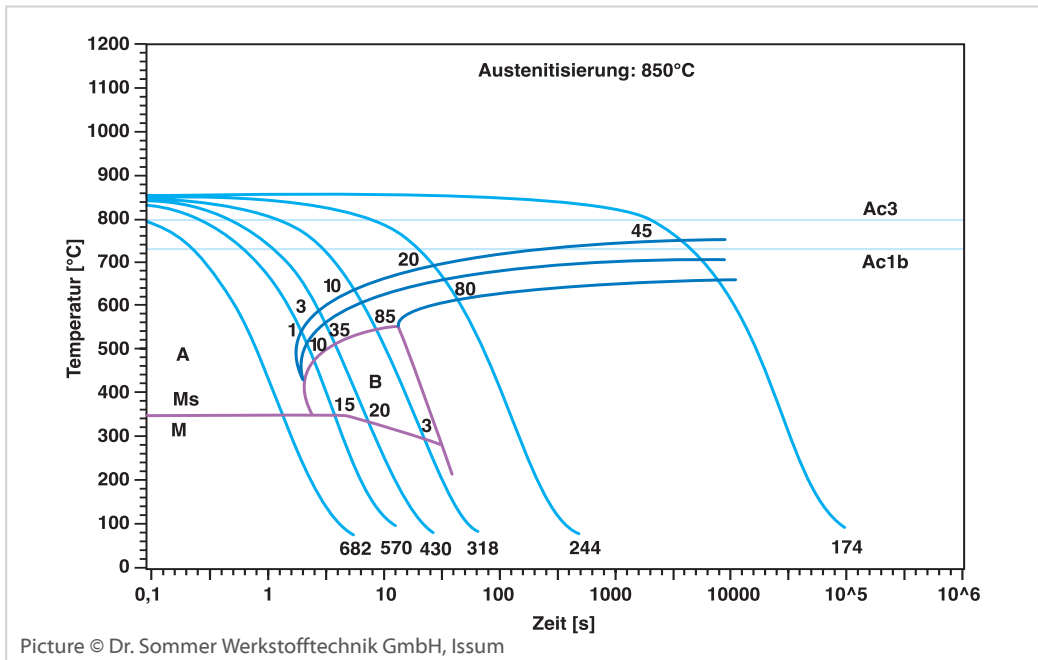

The inductive heating method of the LINSEIS RITA Dilatometer allows the usage of very fast heating and cooling rates. One steel sample like the stainless steel in the application example can be used to determine the CCT, TTT and CHT within only a few measurements. The diagram shows different cooling speeds from 850°C to Room Temperature. The same sample was cooled down several times and for each resulting curve the phase transition points were determined. The purple and dark blue curves are connecting those phase transition points and give the CCT curve for the sample.

(please look at our separate L78 Quenching & Deformation Dilatometer brochure)

Picture © Dr. Sommer Werkstofftechnik GmbH, Issum

**LINSEIS GmbH Germany**

Vielitzerstr. 43

95100 Selb

Tel.: (+49) 9287-880 - 0

Fax: (+49) 9287-70488

E-mail: info@linseis.de

**LINSEIS Inc. USA**

109 North Gold Drive

Robbinsville, NJ 08691

Tel.: (+1) 609 223 2070

Fax: (+1) 609 223 2074

E-mail: info@linseis.com

**LINSEIS China**

Kaige Scientific Park 2653 Hunan Road

201315 Shanghai

Tel.: +86-21-6190 1202

Fax.: +86-21-6806 3576

**LINSEIS France**

Bureaux Paris

52 Boulevard Sébastopol

75003 Paris

Tel.: (+33) 1 73.02.82.72

Bureaux Lyon

2A Chemin des Eglantines

69580 Sathonay Village

(+33) 04.26.68.70.85

**LINSEIS Poland**

ul. Dabrowskiego 1

05-800 Pruszków

Tel.: (+48) 692-773-795

www.linseis.com**Products:** DIL, TG, STA, DSC, HDSC, DTA, TMA, MS/FTIR, In-Situ EGA, Laser Flash, Seebeck Effect, Thin Film Analyzer, Hall-Effect**Services:** Service Lab, Calibration Service

03/18

LINSEIS